


# Scan with confidence

Venue Family point of care ultrasound systems now feature Caption Guidance

## Guiding you to acquire quality cardiac ultrasound images

Now there's a built-in tool to help on the Venue™ Family of point of care ultrasound systems. Thanks to Caption Guidance™ AI-driven software, even new ultrasound users can capture cardiac images successfully. Real-time, turn-by-turn on-screen guidance prompts your probe movements to help capture a diagnostic-quality image.

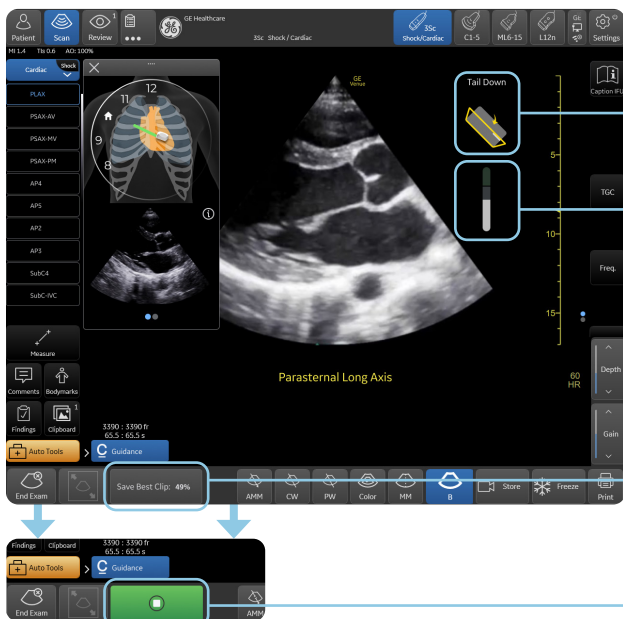

### 1 Prescriptive Guidance

See real-time, turn-by-turn guidance that prompts your movements to properly position your probe for the view you want to capture.

### 2 Quality Meter

Watch the meter rise as the image improves and gets closer to diagnostic-quality.

### 3 Save Best Clip

At the touch of a button, Save Best Clip will select the best available image loop from the continuous recording for each view.

### 4 AutoCapture

For cases where you can reach the quality threshold, AutoCapture will automatically save your clip, hands-free. No need to press record to capture your images.



### Scan your way

Customizable scanning protocols allow you the flexibility to select key cardiac views – up to 10. While scanning, you can easily skip to the view you desire.

### Focused on cardiac

Compatible with the 3Sc-RS phased-array probe and the M5Sc-RS\* broad spectrum sector probe designed for high quality cardiac imaging.



## Helps you before, during and after the scan



#### Before the exam:

Preliminary probe placement information shows you where to place and how to orient the ultrasound probe.

As you begin the exam, a handy reference image shows you the image you are trying to acquire.



#### During the scan:

Follow the prescriptive guidance provided on the screen. Move the probe slowly – and watch the Quality Meter rise as you get closer to an image of diagnostic quality.



#### After scanning:

Easily review your recorded clips prior to ending your exam to make sure you have all the views you need.



\*Only available on Venue system.

Caption Guidance is available in USA only.

Products mentioned in the material may be subject to government regulations and may not be available in all countries. Shipment and effective sale can only occur after approval from the regulator. Please check with local GE HealthCare representative for details.

©2023 GE HealthCare. Venue and Caption Guidance are trademarks of GE HealthCare. GE is a trademark of General Electric Company used under trademark license.

September 2023 JB25668XX